



The POWER to Restore Flow

CLEANERXT™ AND CLEANER15™ Rotational Thrombectomy Systems


ARGON
MEDICAL DEVICES

The Argon Medical Devices family of thrombectomy products





Combining **Torque** with **Trackability**

- **NEW** - Dual power cells increase duration and run time, maintaining torque through tortuosity
- **NEW** - Improved drive shaft design increases efficiency of power transfer and torque
- 6F low profile design enables excellent steerability through vessel tortuosity

The **POWER** to Restore Flow

- Increased sinusoidal wire amplitude – from 9 mm to 15 mm
- Increased wire diameter – from .035” to .044”
- Dual power cells maintain torque through tortuosity


Power Comparison	Power Cells	Voltage	RPM
Cleaner XT™	2	9 Volts	4000
Cleaner 15™	2	9 Volts	4000
Arrow Trerotola™ PTD®	1	6 Volts	3000

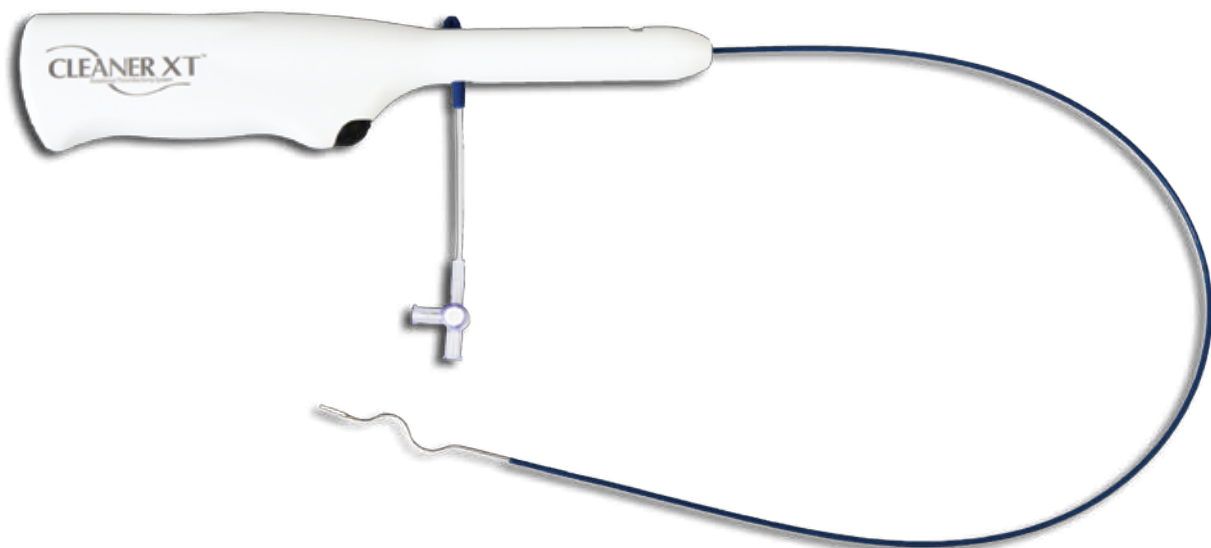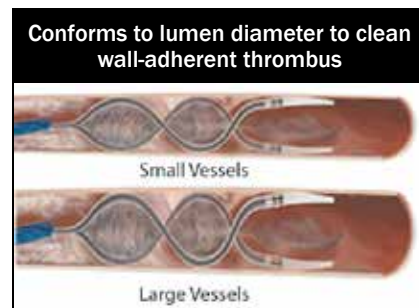

CLEANER XT and CLEANER 15 are indicated for mechanical declotting of native vessel dialysis fistulae and synthetic dialysis access grafts


Innovative Sinusoidal Wire Design

- Radiopaque S-wire conforms to varying lumen diameters, actively cleaning wall-adherent thrombus
- Fluid vortex gently macerates thrombus while reducing the risk of endothelial damage¹
- S-wire is capable of cleaning around opposing sheaths reducing procedure time and sheath manipulation

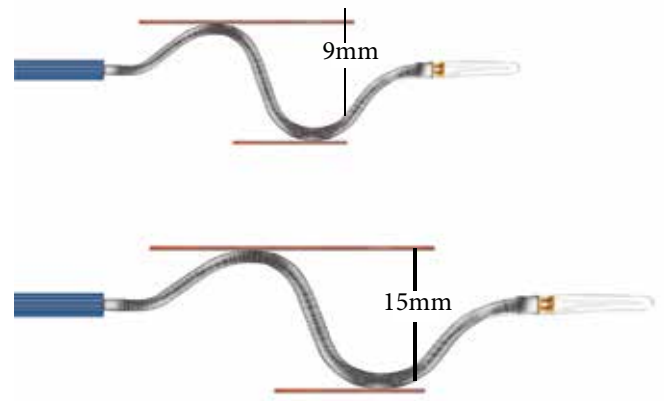

Atraumatic Mechanical Thrombectomy


- Safe and effective in both AV fistulae and grafts
- Guidewire-like design provides excellent steerability
- Available in 65 cm and 135 cm catheter lengths
- 3-way sideport and distal side hole on catheter allows for infusion of fluids and contrast media

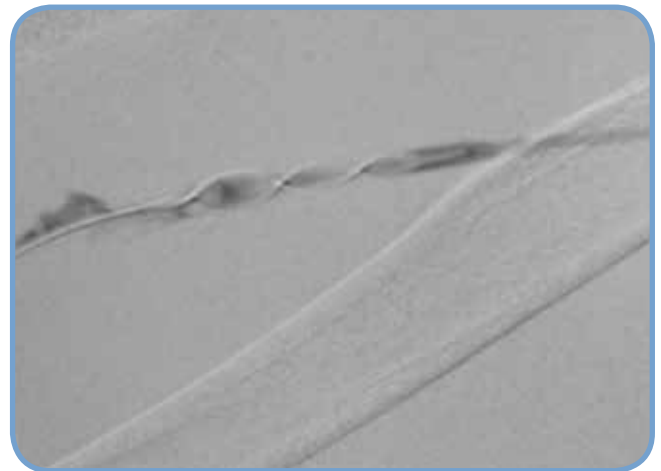

Used with permission of Mark Goldberg, MD

Ergonomic One-Piece Design

- Procedure ready – no assembly required
- Intuitive to operate
- Quiet operation allows for patient comfort during procedure
- Completely disposable – no capital equipment required

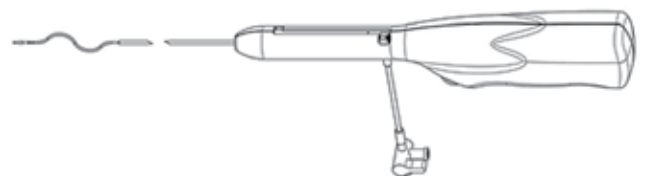

For more information visit www.argonmedical.eu.com



CLEANER XT PRODUCT INFORMATION

Catalogue No.	Catheter OD	Length	Sinusoidal Wire Amplitude
700009XT	6F	65cm	9mm
700109XT	6F	135cm	9mm



CLEANER 15 PRODUCT INFORMATION

Catalogue No.	Catheter OD	Length	Sinusoidal Wire Amplitude
700015	7F	65cm	15mm
700115	7F	135cm	15mm

EMEA Customer Service

Argon Medical Devices Sarl
Espace Performance
Alphasis Bâtiment J
35760 St Gregoire
France
Tel: +33 (0) 2 99 68 94 95
Fax: +33 (0) 2 99 23 87 21
info-europe@argonmedical.com

Direct Sales Denmark

Argon Medical Devices Danmark A/S
Knud Bro Alle 4
DK-3660 Stenløse
Denmark
Tel: +45 47 10 10 48
Fax: +45 70 40 45 10
customer@argonmedical.com

Direct Sales Finland

Argon Medical Devices Danmark A/S
Arenavägen 41 plan 12
SE-121 77 Stockholm
Sweden
Tel: 080030111
Fax: 080030112
customer@argonmedical.com

Direct Sales France

Argon Medical Devices Sarl
Espace Performance
Alphasis Bâtiment J
35760 St Gregoire
France
Tel: +33 (0) 2 99 68 99 24
Fax: +33 (0) 2 99 23 87 21
service-france@argonmedical.com

Direct Sales Germany

Argon Medical Devices Netherlands B.V.
Zweigniederlassung Deutschland
Steinpilzweg 26
91058 Erlangen
Deutschland
Tel: +49 (0)9131 9085969
Fax: +49 (0)9131 9081964
orders.gsa@argonmedical.com

Direct Sales Norway

Argon Medical Devices Danmark A/S
Knud Bro Alle 4
DK-3660 Stenløse
Denmark
Tel: +47 80010046
Fax: +45 80187692
customer@argonmedical.com

Direct Sales Sweden

Argon Medical Devices Danmark A/S
Arenavägen 41 plan 12
SE-121 77 Stockholm
Sweden
Tel: +46 8 50 56 88 70
Fax: +46 8 50 56 88 69
customer@argonmedical.com



www.argonmedical.eu.com

© 2015 Argon Medical Devices, Inc., All Rights Reserved.

No reproduction of any portion of this work shall be made without prior permission of the author.

96-6300-43/0615A